

FAM-VAD-FMK

 Catalog number: 13470, 13484
 Unit size: 25 Tests, 100 Tests

Component	Storage	Amount (Cat No. 13470)	Amount (Cat No. 13484)
FAM-VAD-FMK	Freeze (< -15 °C), Minimize light exposure	25 Tests	100 Tests

OVERVIEW

FAM-VAD is a green fluorescent cell-permeable polycaspase inhibitor to target caspases 1, 2, 3, 6, 8, 9, or 10. Once inside the cell, the inhibitor binds covalently to the active caspase, which produces a green fluorescence. When added to a population of cells, the FAM-VAD-FMK probe enters each cell and covalently binds to a reactive cysteine residue that resides on the large subunit of the active caspase heterodimer, thereby inhibiting further enzymatic activity. The bound labeled reagent is retained within the cell, while any unbound reagent will diffuse out of the cell and is washed away. The green fluorescent signal is a direct measure of the amount of active caspase present in the cell at the time the reagent was added. Cells that contain the bound labeled reagent can be analyzed by 96-well plate-based fluorometry, fluorescence microscopy, or flow cytometry. It works in diverse cell lines such as human, rodent and drosophila.

AT A GLANCE

Important It is important to store at <-15 °C and should be stored in cool, dark place.

Following protocol only provides a guideline, and should be modified according to your specific needs.

It can be used within 12 months from the date of receipt.

SAMPLE EXPERIMENTAL PROTOCOL
General Solution Caspase Assays Using AMC, AFC, pNA, R110 and ProRed Substrates

1. Prepare a 10 mM stock solution in DMSO.
2. Prepare a 2X caspase substrate (50 μM) assay solution as the following: 50 μL substrate stock solution, 100 μL DTT (1M), 400 μL EDTA (100 mM), 10 mL Tris Buffer (20 mM), pH =7.4.
3. Mix equal volume of the caspase standards or samples with 2X caspase substrate assay solution, and incubate the solutions at room temperature for at least 1 hour.
4. Monitor the fluorescence using a fluorescence microplate reader, or absorbance using an absorbance microplate reader.

Cell Caspase Assays Using Cell-Permeable FMK Caspase Probes

1. Prepare a 2-5 mM stock solution in DMSO.
2. Treat cells as desired.
3. Prepare a 2X permeable caspase substrate (20 μM) assay solution by diluting the DMSO stock solution in Hanks with 20 mM Hepes buffer (HHBS).
4. Mix equal volume of the treated cells with 2X caspase substrate assay solution, and incubate the cells in a 37 °C, 5% CO₂ incubator for at least 1 hour.
5. Wash the cells with HHBS for at least once.
6. Monitor the fluorescence intensity by a flow cytometer, a fluorescence microscope or a fluorescence microplate reader.

Cell Caspase Assays Using Cell-Permeable FMK Caspase Probes (For #13470-13476 only)

1. Prepare a 250X stock solution by adding 50 μL DMSO into the vial.
Note For Cat# 13484, prepare a 250X stock solution by adding 200 μL DMSO into the vial.
2. Treat cells as desired.
3. Add 250 X DMSO stock solution into the cell solution at a 1:250 ratio (such as 2 μL to 500 μL cells), and incubate the cells in a 37 °C, 5% CO₂ incubator for 1 hour.
4. Wash the cells with HHBS for at least once.
5. Monitor the fluorescence intensity by flow cytometer, fluorescence microscopy or fluorescent microplate reader.

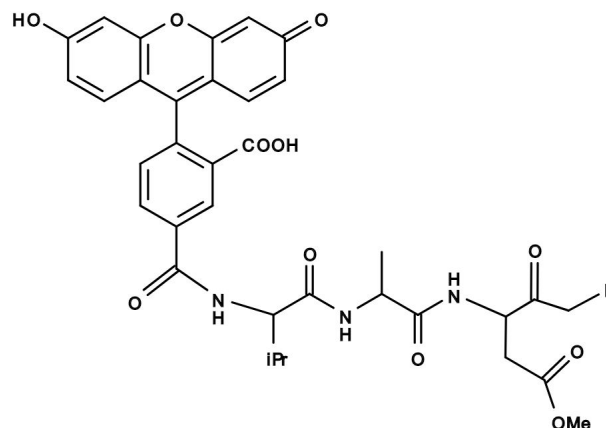
EXAMPLE DATA ANALYSIS AND FIGURES


Figure 1. Chemical structure for FAM-VAD-FMK

DISCLAIMER

AAT Bioquest provides high-quality reagents and materials for research use only. For proper handling of potentially hazardous chemicals, please consult the Safety Data Sheet (SDS) provided for the product. Chemical analysis and/or reverse engineering of any kit or its components is strictly prohibited without written permission from AAT Bioquest. Please call 408-733-1055 or email info@aatbio.com if you have any questions.