

TAMRA-cGMP PDE V substrate *Red Fluorescence*

 Catalog number: 13605
 Unit size: 0.5 umol

Component	Storage	Amount
TAMRA-cGMP PDE V substrate *Red Fluorescence*	Freeze (< -15 °C), Minimize light exposure	1 vial (0.5 umol)

OVERVIEW

This red cGMP derivative is a specific substrate for phosphodiesterase (PDE) V. It can be used for assaying PDE V activities or screening PDE V inhibitors in combination with anti-cGMP antibody in a FRET readout or FP format. PDE is a group of enzymes that degrade the second messenger molecules: cyclic nucleotides cAMP and cGMP. They regulate the localization, duration, and amplitude of cyclic nucleotide signaling within subcellular domains. PDEs are therefore important regulators of signal transduction mediated by these second messenger molecules. PDE enzymes are often targets for pharmacological inhibition due to their unique tissue distribution, structural and functional properties. Inhibitors of PDE can prolong or enhance the effects of physiological processes mediated by cAMP or cGMP by inhibition of their degradation by PDE. PDE inhibitors have been identified as new potential therapeutics in areas such as pulmonary arterial hypertension, coronary heart disease, dementia, depression and schizophrenia. For example, Sildenafil (Viagra) is an inhibitor of cGMP-specific PDE V, which enhances the vasodilatory effects of cGMP in the corpus cavernosum, and is used to treat erectile dysfunction.

AT A GLANCE

Important Following protocol only provides a guideline, and should be modified according to your specific needs.

KEY PARAMETERS
Fluorescence microplate reader

Excitation	540 nm
Emission	590 nm
Cutoff	570 nm
Recommended plate	Solid black

PREPARATION OF STOCK SOLUTIONS

Unless otherwise noted, all unused stock solutions should be divided into single-use aliquots and stored at -20 °C after preparation. Avoid repeated freeze-thaw cycles.

TAMRA-cGMP PDE V stock solution (1 mM)

Make a 1 mM stock solution by adding 500 µL of DMSO into the vial of 0.5 umol TAMRA-cGMP PDE V substrate.

PREPARATION OF WORKING SOLUTION
TAMRA-Cyclic-3', 5'-GMP PDE V substrate assay solution (2X)

Make 2X TAMRA-Cyclic-3', 5'-GMP PDE V substrate assay solution by diluting 1 mM TAMRA-Cyclic-3', 5'-GMP PDE V substrate stock solution into your PDE buffer (such as 10 mM Tris-HCl, pH 7.4, 10 mM Mg Cl₂, 1 mM MnCl₂) to make a 200 - 400 nM solution.

Note Make only sufficient quantity needed for the assay.

SAMPLE EXPERIMENTAL PROTOCOL

- Mix equal volume of the PDE V standards or samples with 2X TAMRA-Cyclic-3', 5'-GMP PDE V substrate assay solution, and incubate at room temperature for at least 1 hour.

- Monitor the fluorescence polarization at Ex/Em = 540/590 nm.

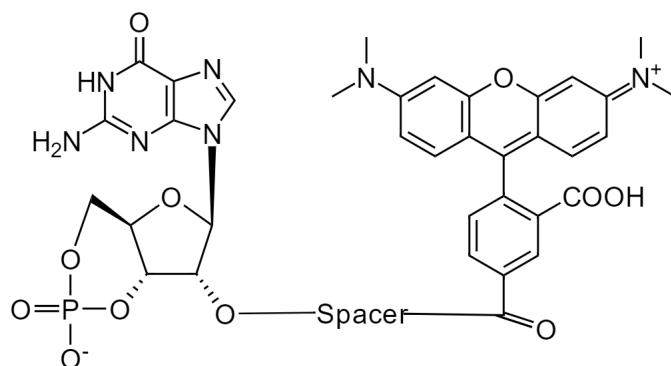
EXAMPLE DATA ANALYSIS AND FIGURES


Figure 1. Chemical structure for TAMRA-cGMP PDE V substrate *Red Fluorescence*

DISCLAIMER

AAT Bioquest provides high-quality reagents and materials for research use only. For proper handling of potentially hazardous chemicals, please consult the Safety Data Sheet (SDS) provided for the product. Chemical analysis and/or reverse engineering of any kit or its components is strictly prohibited without written permission from AAT Bioquest. Please call 408-733-1055 or email info@aatbio.com if you have any questions.