

Mag-520™ AM

 Catalog number: 20406
 Unit size: 10x50 ug

Component	Storage	Amount
Mag-520™ AM	Freeze (< -15 °C), Minimize light exposure	10x50 ug

OVERVIEW

Intracellular magnesium is important for mediating enzymatic reactions, DNA synthesis, hormone secretion, and muscle contraction. A vast majority of the existing magnesium ion indicators are based on tricarboxylate APTRA chelator derived from the popular tetracarboxylate BAPTA calcium ion chelator. They include mag-fura-2, mag-indo-1, mag-fluo-4 and mag-rhod-3. However, all of them have higher affinity for calcium than magnesium although they were designed for detecting magnesium ion. Mag-520™ is the first commercial magnesium indicator that has higher affinity for magnesium than calcium. Its significantly improved selectivity can be used for magnesium signaling applications with fluorescence microscopy or flow cytometry. Mag-520™ AM is cell-permeable with almost identical spectra to those of FITC, making it a convenient fluorescent probe used with FITC filter set that is ubiquitously equipped with almost all fluorescence instruments.

AT A GLANCE

Protocol summary

1. Grow cells as desired
2. Prepare and add Mag-520™ AM working solution to samples
3. Incubate samples at 37 °C for 15 to 45 minutes
4. Monitor the fluorescence intensity using flow cytometer with 530/30 nm filter (FITC channel) or using fluorescence microscopy with FITC filter set

KEY PARAMETERS

Flow cytometer

Excitation 488 nm laser
 Emission 530/30 nm filter
 Instrument specification(s) FITC channel

Fluorescence microscope

Excitation FITC filter set
 Emission FITC filter set
 Recommended plate Black wall/clear bottom

PREPARATION OF STOCK SOLUTIONS

Unless otherwise noted, all unused stock solutions should be divided into single-use aliquots and stored at -20 °C after preparation. Avoid repeated freeze-thaw cycles.

Mag-520™ AM stock solution

Add appropriate amount of DMSO into Mag-520™ AM vial to make 2-5 mM Mag-520™ AM stock solution.

Note Mag-520™ AM working solution should be used promptly.

PREPARATION OF WORKING SOLUTION

Mag-520™ AM working solution

Prepare 5-50 μM Mag-520™ AM working solution in buffer of your choice.

Note Mag-520™ AM working solution should be used promptly.

Note The concentration of the Mag-520™ AM should be optimized for different cell types and conditions.

SAMPLE EXPERIMENTAL PROTOCOL

The following protocol can be used as a guideline and should be optimized according to the needs.

1. Grow cells as desired.
2. Remove the treatment and wash the cells with buffer of your choice such as DPBS.
3. Add Mag-520™ AM working solution and incubate the samples for 15-45 minutes at 37 °C incubator.

Note Optimal time for incubation needs to be determined carefully.

4. Remove the working solution and wash cells with buffer of your choice.
5. Resuspend cells in buffer of your choice.
6. Add stimulant to stimulate the cells and monitor the fluorescence intensity with flow cytometer using 530/30 nm filter (FITC channel) or fluorescence microscope with FITC filter set.

EXAMPLE DATA ANALYSIS AND FIGURES

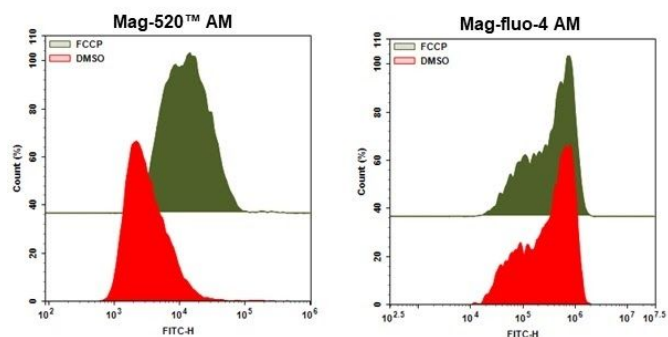


Figure 1. In vitro detection of Mg⁺². HL-60 cells were stained with AM form of Mag-520™ (5 μM) and mag-fluo-4 (5 μM) for 30 minutes and stimulated with FCCP (20 μM) for 30 minutes and response was recorded for end point analysis. Response was recorded using NovoCyte™ 3000 Flow Cytometer with 530/30 nm filter and 488 nm laser.

DISCLAIMER

AAT Bioquest provides high-quality reagents and materials for research use only. For proper handling of potentially hazardous chemicals, please consult the

Safety Data Sheet (SDS) provided for the product. Chemical analysis and/or reverse engineering of any kit or its components is strictly prohibited without written permission from AAT Bioquest. Please call 408-733-1055 or email info@aatbio.com if you have any questions.