

## ReadiUse™ Preactivated APC-iFluor™ 800 Tandem

 Catalog number: 2572  
 Unit size: 1 mg

Component	Storage	Amount
A: ReadiUse™ Preactivated APC-iFluor™ 800 Tandem	Refrigerated (2-8 °C), Minimize light exposure	1 vial (1 mg)
B: Buccutite™ MTA	Freeze (< -15 °C), Minimize light exposure	1 vial (100 µg)
C: Spin Desalting Column		Not Included

### OVERVIEW

Allophycocyanin (APC) is a phycobiliprotein isolated from *Spirulina* sp., a blue-green alga. Like other phycobiliproteins, APC is fluorescent, with an extremely high absorptivity and a high quantum efficiency. Our APC-iFluor™ 800 Tandem is a unique fluorophore with an IR emission developed for spectral flow cytometry. Its emission spectrum does not overlap with any of the common fluorophores, making a unique color for multicolor applications in flow cytometry. ReadiUse™ Preactivated APC-iFluor™ 800 Tandem is an activated APC protein, and can be easily conjugated to antibodies with much higher conjugation yield than the conventional APC tandems.

### SAMPLE EXPERIMENTAL PROTOCOL

#### Preparation of pre-activated antibody with Buccutite™ MTA

1. Reconstitute Buccutite™ MTA in DMSO at ~10 mg/mL.
 

**Note** Please store unused Buccutite™ MTA at -20 °C and could be used up to two freeze and thaw cycles.
2. Prepare target antibody (Ab) in pH = 8.5 - 9.0 buffer at concentration above 1 mg/mL.
3. Add Buccutite™ MTA to Ab solution at the ratio of 8 - 10 µg Buccutite™ MTA/100 µg Ab.
4. Mix well and react at room temperature for 60 minutes, rotating during the reaction.
5. Purify the reaction mixture with desalting column to remove unreacted Buccutite™ MTA and exchange buffer to PBS or buffer of your choice.
6. Collect the Buccutite™ MTA-activated Ab, and estimate the concentration by 70% yield of the original starting amount.

#### Conjugation with pre-activated APC-iFluor™ 800 Tandem


1. Reconstitute pre-activated APC-iFluor™ 800 Tandem in 100 µL ddH<sub>2</sub>O to 10 mg/mL.
 

**Note** Reconstituted pre-activated APC-iFluor™ 800 Tandem could be stored at 4 °C for one month, kept from light.
2. Add APC-iFluor™ 800 Tandem directly to MTA-activated target Ab solution at the ratio of 130 µg APC-iFluor™ 800 Tandem/100 µg MTA-activated Ab.
3. Rotate the mixture for 60 minutes at room temperature.
4. The Ab/APC-iFluor™ 800 Tandem conjugates are now ready to use.

**Note** The antibody conjugate should be stored at > 0.5 mg/mL in the presence of a carrier protein (e.g., 0.1% bovine serum albumin) and 0.02% - 0.05% sodium azide. The Ab/APC-iFluor™ 800 Tandem conjugates solution could be stored at 4 °C for up to two months, and kept from light.

5. (Optional) Ab/APC-iFluor™ 800 Tandem conjugates could be further purified through size exclusion chromatography to get best performance.

### EXAMPLE DATA ANALYSIS AND FIGURES

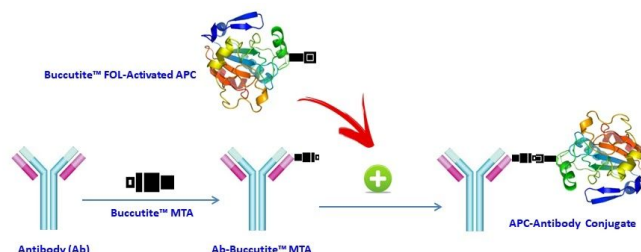

**Figure 1.** Our preactivated APC was premodified with our Buccutite™ FOL (provided). Your antibody (or other proteins) is modified with our Buccutite™ MTA (provided as free sample) to give MTA-modified protein (such as antibody). The MTA-modified protein readily reacts with FOL-modified APC (provided) to give the desired APC-antibody conjugate in much higher yield than the SMCC chemistry. In addition our preactivated APC reacts with MTA-modified biopolymers at much lower concentrations than the SMCC chemistry.

### DISCLAIMER

AAT Bioquest provides high-quality reagents and materials for research use only. For proper handling of potentially hazardous chemicals, please consult the Safety Data Sheet (SDS) provided for the product. Chemical analysis and/or reverse engineering of any kit or its components is strictly prohibited without written permission from AAT Bioquest. Please call 408-733-1055 or email [info@aatbio.com](mailto:info@aatbio.com) if you have any questions.