

Cell Meter™ Cellular Senescence Activity Assay Kit

Catalog number: 23005
Unit size: 100 Tests

Component	Storage	Amount
Component A: Xite™ beta-D-galactopyranoside	Freeze (<-15 °C), Minimize light exposure	1 vial
Component B: Assay Buffer	Freeze (<-15 °C)	1 bottle (20 mL)
Component C: DMSO	Freeze (<-15 °C)	1 vial (100 uL)

OVERVIEW

Cellular Senescence is an irreversible growth arrest triggered in order to prevent growth in DNA damaged cells. Senescence-associated beta-galactosidase (SA-beta-gal) is highly overexpressed in senescent cells and it has been widely used as a senescence marker. X-gal staining, a colorimetric method is widely available and used to detect SA-beta-gal in senescent cells. The color method has some limitations such as requirement of fixation of samples due to the low cell permeability of X-gal, longer staining time and low sensitivity. Cell Meter™ Cellular Senescence Activity Assay Kit uses Xite™ beta-D-galactopyranoside, a fluorogenic beta-Gal substrate that readily enters into live cells, and gets cleaved by SA-β-gal inside cells, generating strong green fluorescence. Unlike cell-impermeable X-Gal substrate, it has excellent cell permeability. Cell Meter™ Cellular Senescence Activity Assay Kit enables users to detect the senescence with higher sensitivity with robust performance. The Xite product is well retained inside the cells, producing a stable signal for fluorescence imaging and flow cytometry analysis.

AT A GLANCE

Protocol summary

1. Treat samples as desired
2. Prepare and add Xite™ beta-D-galactopyranoside working solution to samples
3. Incubate samples at 37 °C for 15 to 45 minutes
4. Monitor the fluorescence intensity using flow cytometer with 530/30 nm filter (FITC channel)

Important Thaw each kit component at room temperature before starting the experiment.

KEY PARAMETERS

Instrument: Flow cytometer
Excitation: 488 nm laser
Emission: 530/30 nm filter
Instrument specification(s): FITC channel

Instrument: Fluorescence microscope
Excitation: FITC filter set
Emission: FITC filter set
Recommended plate: Black wall/clear bottom

PREPARATION OF STOCK SOLUTIONS

Unless otherwise noted, all unused stock solutions should be divided into single-use aliquots and stored at -20 °C after preparation. Avoid repeated freeze-thaw cycles.

Xite™ beta-D-galactopyranoside stock solution (100X):

Add 100 uL DMSO (Component C) into Xite™ beta-D-galactopyranoside (Component A) and mix well.

Note Store the unused Xite™ beta-D-galactopyranoside stock solution at -20 °C in single use aliquots.

PREPARATION OF WORKING SOLUTION

Xite™ beta-D-galactopyranoside working solution (1X):

Dilute 10 uL of Xite™ beta-D-galactopyranoside stock solution (100X) with 1 mL of Assay Buffer to make Xite™ beta-D-galactopyranoside working solution (1X).

Note Xite™ beta-D-galactopyranoside working solution should be used promptly.

SAMPLE EXPERIMENTAL PROTOCOL

The following protocol can be used as a guideline and should be optimized according to the needs.

1. Treat your samples as desired.
2. Wash the cells with buffer of your choice such as DPBS.
3. Add 100 uL Xite™ beta-D-galactopyranoside working solution for 15-45 minutes and incubate the samples at 37°C incubator.

Note Optimal time for incubation needs to be determined carefully.

4. Remove the working solution and wash cells with buffer of your choice.
5. Resuspend the cells in the Assay Buffer (Component B) and monitor the fluorescence intensity with flow cytometer using 530/30 nm filter (FITC channel) or fluorescence microscope with FITC filter set.

EXAMPLE DATA ANALYSIS AND FIGURES

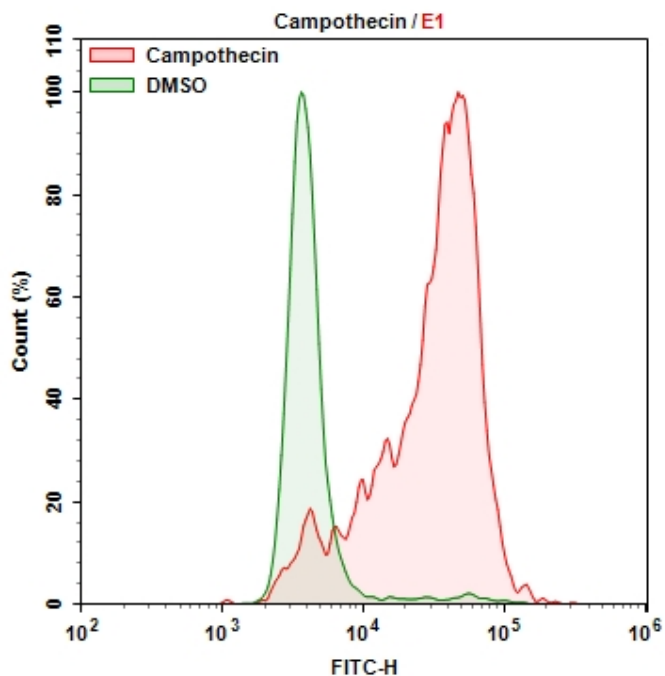


Figure 1. Cellular senescence was measured with Cell Meter™ Cellular Senescence Activity Assay Kit using a NovoCyte Flow Cytometer (ACEA Biosciences). HL-60 cells were incubated with Camptothecin for 6 hours to induce senescence and stained with Xite™ beta-D-galactopyranoside for 30 mins at 37°C. The signal was acquired using FITC channel in ACEA NovoCyte flow cytometer.

DISCLAIMER

AAT Bioquest provides high-quality reagents and materials for research use only. For proper handling of potentially hazardous chemicals, please consult the Safety Data Sheet (SDS) provided for the product. Chemical analysis and/or reverse engineering of any kit or its components is strictly prohibited without written permission from AAT Bioquest. Please call 408-733-1055 or email info@aatbio.com if you have any questions.